

7in1 Hydra facial machine



User manual

USER MANUAL MULTI FUNCTION SKIN CARE SYSTEM

IMPORTANT: before starting the equipment correctly, please follow the installation instructions. Any damages caused by incorrect operations will void the warranty.

SAFETY INFORMATION

Thank you for choosing our products! This user instruction is specially designed for the safety of you and others. Please read the instructions carefully before installing or using the machine and follow them when using.

1.1 Electrical Safety

- Before installing the machine, make use that your local electric network complies with the voltage and frequency indicated on the tag on the rear panel. Otherwise, the machine may be damaged in case of excessive current or voltage rise.
- Never touch the machine with damp or wet hands or feet.
- It is dangerous to change or try to modify the function of this machine in any way.
- Do not attempt to operate this unit if the power cord has been frayed or broken
- Do not spill water or other liquids into or on to your unit.
- During long periods of non-use, disconnect the unit's main power.
-


1.2 Safety Rules When Installing

- Place the machine on a stable and secure base. It should not be subjected to any jolts or vibrations.
- Never expose the unit to too high ($>40^{\circ}\text{C}$) nor too low temperature ($<5^{\circ}\text{C}$).
- NEVER EXPOSE THE UNIT TO RAIN OR ANY OTHER HUMIDITY.
- Do not expose the unit to direct sunshine or install the unit close to a source of heat (central heating, radiator...); otherwise the screen may be damaged.
- Do not cover or block the ventilator behind the unit.
-



1.3 Safety Rules When in Use.

- Well clean areas to treat before treatment.
- Never touch the machine with damp or wet hands or feet.

I. Instrument display

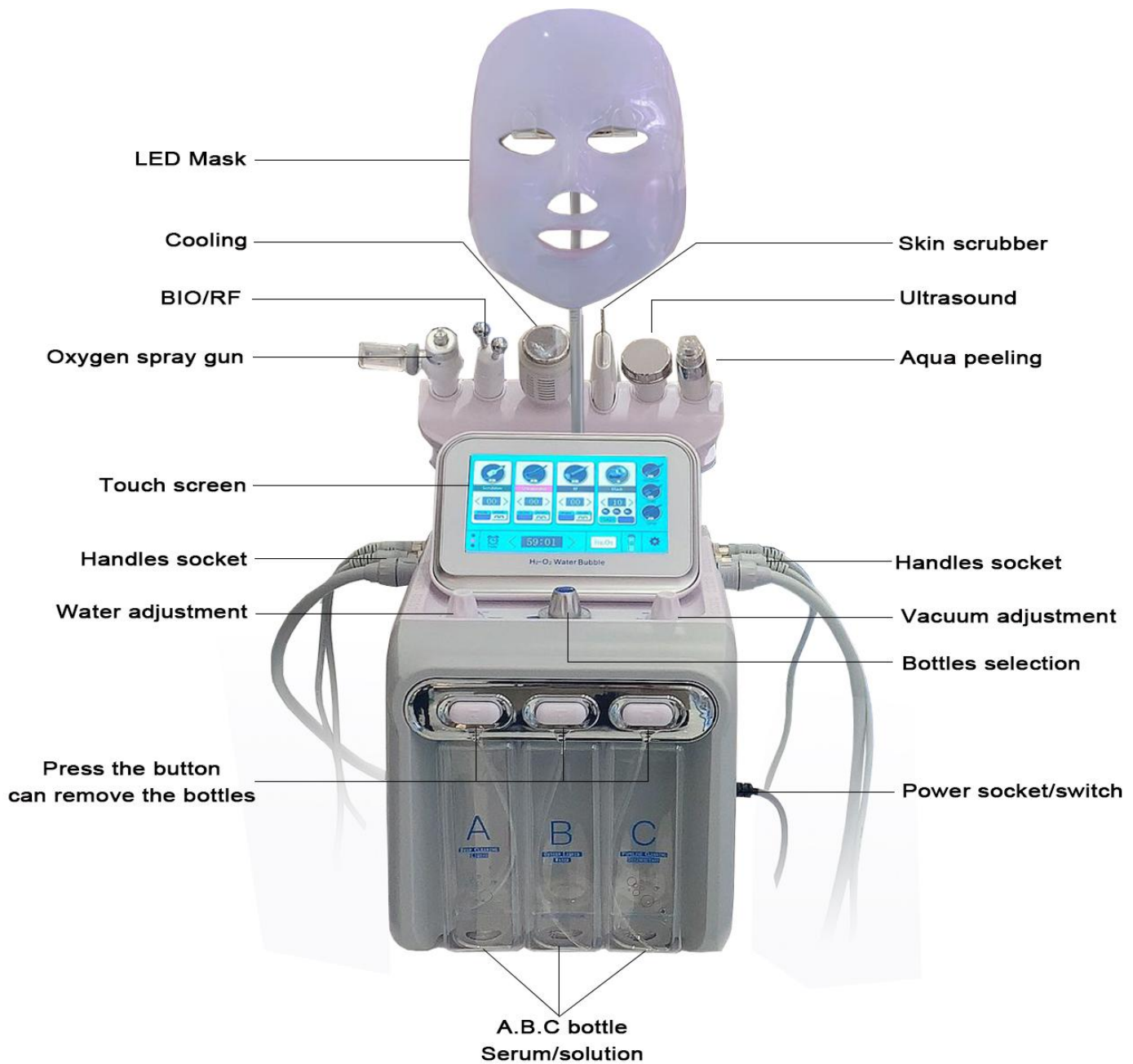

A.B.C bottle addable product solution(Such us:cleans skin solution, soothes water, acne skin water and pure water etc)

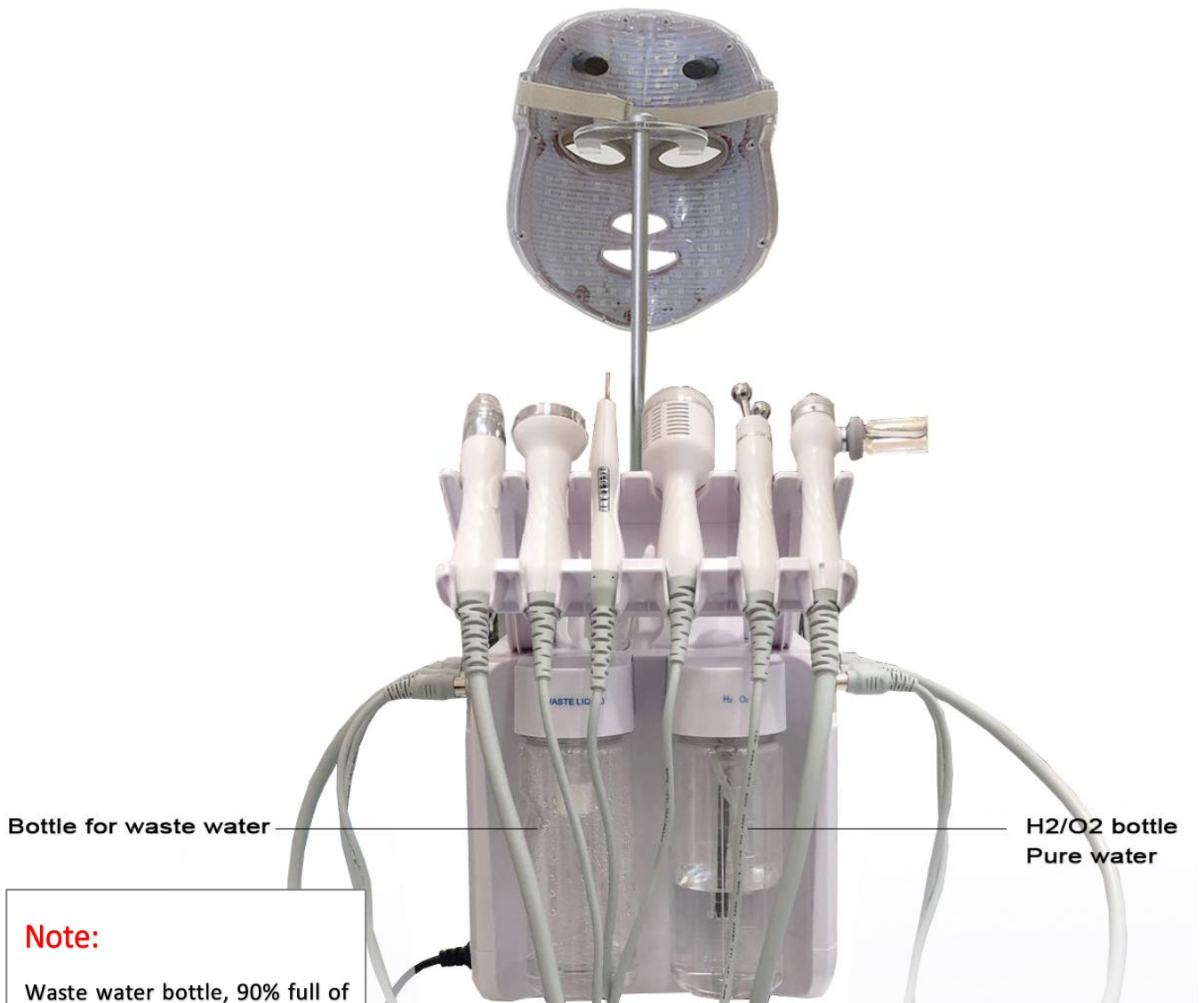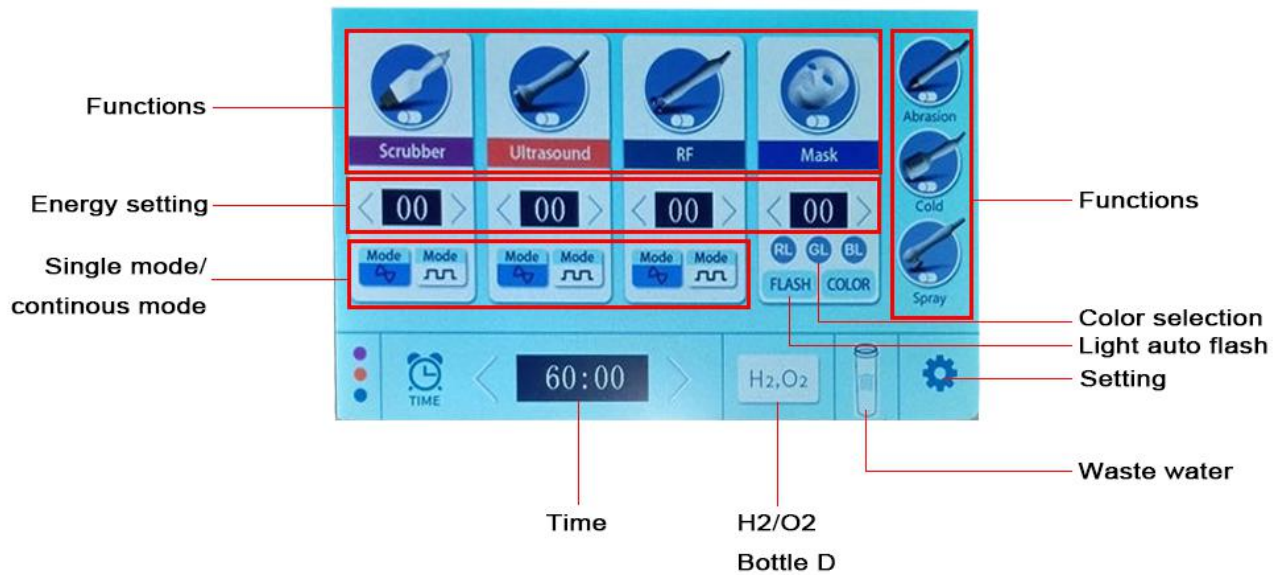

Note:


Waste water bottle, 90% full of water must be dumped, otherwise it will damage the instrument

	Sprayer: High pressure oxygen injection makes skin absorb moisture more easily		Repair skin, improving pigmentation, delaying aging, etc
	Cooling: Shrink pores Firming skin		Ultrasonic: promote Micro circulation, skin rejuvenating
	BO lifting: raise the overall tightening of facial skin		Aqua peeling: Deep cleansing and refreshing
	Scrubber: facial cleansing, peeling, pore cleansing		

B. Instrument operation

1. Install all handles to related handle socket. Input solution inside A.B.C bottles. Input pure water inside D bottle. Connect power cord to the power socket. Switch on.

Aqua peeling, select the solution bottle, set the vacuum and water. Press Start on the system or press the button on the handle.

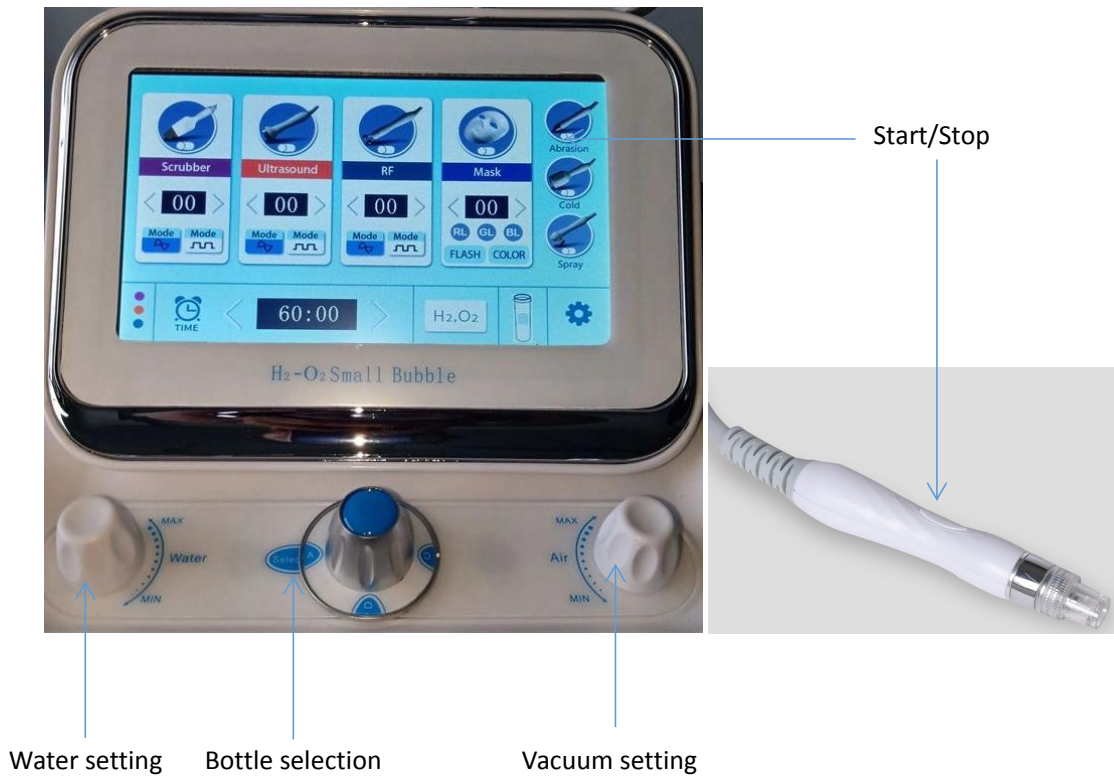

After treatment, select bottle D, turn off H₂/O₂ function. use finger to block the handle tips for 2minutes for cleaning handle pipe

Scrubber, Apply solution on the skin, Set intensity 5, Start on the system, press metal on the handle to remove dead skin, black head.etc



Oxygen spray, Input toner or acid inside the bottle of spray gun. Press Start on the system or press the button on the handle.



After treatment, input warm water inside the bottle of spray gun, and spray all water out for handle clean.


BIO/RF Skin lifting, apply cream or essence on the skin. Start intensity from 3, Press Start on the system or press the button on the handle to start treatment. Recommend 8-10mins



Ultrasound, apply cream or essence on the skin. Start intensity from 3, Press Start on the system or press the button on the handle to start treatment. Recommend 8-10mins



LED mask, clean the skin. Apply paper mask on the face, set mask intensity 5m cover LED mask on the face for 15minutes.

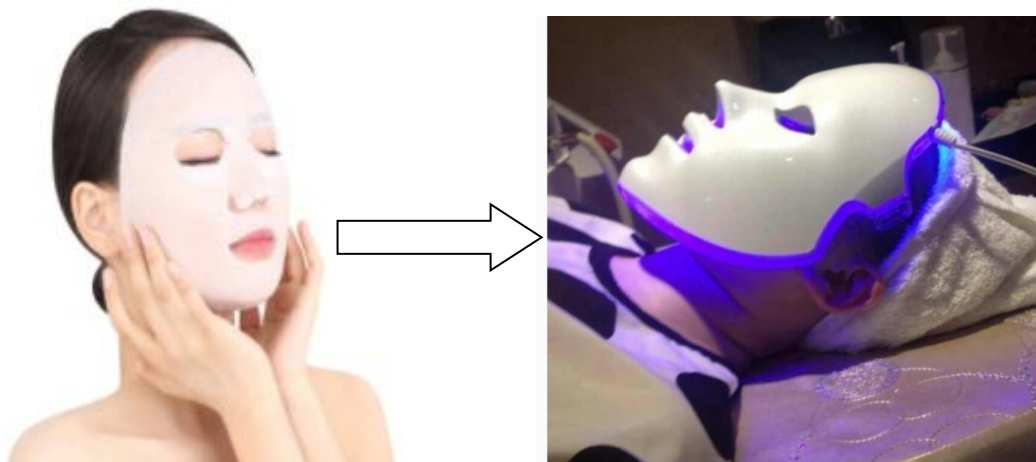

Cooling, After LED mask treatment, use cold hammer to help the skin Soothing, repairing, calming. Recommend 8-10 minutes



Tips:

Hydrogen-oxygen bubble mode: Pour the nutrient solution into A solution bottle or introduce purified water into B solution bottle; start up the ABSORB icon; adjust the working hours, water output, vacuum suction and hot air bubble gear; perform functional tests on the hand to ensure that the suction and water output and temperature are properly re-operated for the customer. Apply the vertical probe of the skin pen to the skin and adsorb.

Ice hammer (COLD), lighting icon can work; high molecular water oxygen (SPRY) after lighting icon can choose low, medium, high gear, ampoule bottles need to use water agent products (e.g. essence, saline, etc.).

Ion rollers (BIO), plasma bipolar desensitization (STRENGTH), etc. Model Light up the icon and choose the working time, mode, energy (from low to high) depending on customer acceptance capacity).

Nano-plastic introduction (DOSE) requires adjustment of dose and energy intensity (from low to high) depending on customer acceptance capacity).

1. Points for attention

Note:

1.1: Before use, please fill the solution bottle with proper pure water or nutrient solution, do not use tap water;

1.2: Pay attention to the level of the waste water bottle when it reaches the upper limit will alarm, please pour it out in time to prevent damage to the instrument.

1.3: Skin cuticle thin or zygomatic red blood, careful use.

1.4: People who have had surgery or deep exfoliation on their faces who haven't reached recovery are careful.

1.5: after operation need to strengthen water supply, a week of continuous application of water surface membrane, do not use efficacy type of water supply products.

2. End of operation

2.1 START/STOP key, then Closed Total power switch. After the operation is finished, gently wipe the control handle with a soft cloth and place it back on the control handle holder.

3. Maintenance and repair

3.1 Maintenance

3.1.1 Routine inspection

The management must inspect the instruments and components.

(a) whether there is a reliable connection between the power outlet and the

ground.

(b) Whether the power line is deformed, broken or broken: by the naked eye, whether the wire is deformed or deformed, if any, it may be caused by leakage of electricity And caused the fire. please replace the new power cord immediately.

c) whether the control handle is abnormal, broken line: through the naked eye inspection, whether the output line has deteriorated, deformation, if any, please replace immediately.

d) Please check that the instrument is correct before use and stop using it.

3.1.2 Cleaning

When cleaning the instrument, the dust must be wiped with a dry soft cloth, and the stolen goods are wiped with a soft wet cloth. Before cleaning this instrument, the power switch must be turned off

3.1.3 The use of the instrument for a long period of time (more than one year) should be carried out by professional maintenance personnel.

3.1.4 When the instrument is not in use, it should be placed in a ventilated and dry place, covered with a dust shield.

3.2 Maintenance

3.2.1 The warranty card must be filled in and sent to the manufacturer or agent after purchase of the instrument to have the warranty effect; the warranty is free of charge for one year from the date of issuing the bill for the quality of the products of our company (Note: with the purchase of the instrument note and the warranty card, the warranty is limited to the instrument mainframe, excluding the spare parts and other consumable items which need to be replaced regularly).

3.2.2 After the warranty period, please contact the manufacturer directly to provide the service guarantee.

3.2.3 The maintenance of this machine must contact our customer service department or dealer, do not repair without authorization, self-opened, the company does not guarantee.

4. Common Troubles hooting Methods

4.1 The instrument cannot work properly, do not easily disassemble and assemble, should be timely with the company's customer service department or dealers, after finding out the cause of failure.

4.2 Before requesting maintenance services, please click the table below to see if you can handle the problem yourself.

5. Responsibility of the manufacturer

5.1 Our company implements life-long maintenance and technical service for our products, and the manufacturing quality problems that occur during the stipulated quality assurance period will be repaired for users free of charge.

5.2 The Company shall not be liable for property damage arising from improper use by the

Customer:

(a) if the purchaser or operator does not follow the instructions and does not follow the safety

precautions or warning instructions;

(b) in violation of the regulations on the transportation, installation, use, maintenance and storage

of the product resulting in damage to the product, person and property.

11. Replacement of accessories

Accessories must be purchased from our company when they need to be replaced.

6. Storage and transportation

6.1 Storage: the instrument should be stored in the original packing box, in a well-ventilated

clean room, the packing box should be high cushion, the ambient temperature is -10°C ~40°C, No more than 93% relative humidity, no harmful gases, flammable, explosive substances and

corrosive gases.

6.2 Transport: Instruments under packing conditions are suitable for road, rail, air and water

transportation and should prevent violent vibration and impulse during handling and

transportation

Strike, and shall not be dampened, and shall not be mixed with flammable, corrosive substances, mixed transport, specific requirements according to the order contract.

7. Duration of use

7.1 The warranty for this instrument is 1 year. (Note: Except for man-made damage)

7.2 The warranty period for accessories is 3 months. (Note: Except for man-made damage)

7.3 The continued use of the instrument beyond its service life shall be subject to the approval of

the Engineer's examination and approval, and the performance of the instrument shall be periodically tested by the Engineer in subsequent use.

7.4 End-of-life instruments and accessories beyond their service life shall be disposed of in accordance with local regulations